

Alta™ Multipurpose Device

AN ALL IN ONE WORKFLOW EFFICIENCY SOLUTION

Alta is an all in one solution that reduces setup times and quickly adapts to a wide range of setups, including SRS, Intracranial, Head & Neck, Craniospinal, Whole Brain, Breast, Pelvis, Abdomen, Lung, Thorax, and SBRT.

- Combines functionality of multiple devices
- Dependable and reproducible
- Fast and efficient setups
- Can be extended beyond end of couchtop for improved head and neck dosimetry
- Lightweight-5.4 kg (12 lbs)
- 61 cm (24 in) W x 199 cm (78.5 in) L
- MR Safe, CT, and treatment compatible



SRS



PELVIS



INTRACRANIAL, HEAD & NECK



ABDOMEN



CRANIOSPINAL



LUNG



WHOLE BRAIN




THORAX



BREAST



SBRT

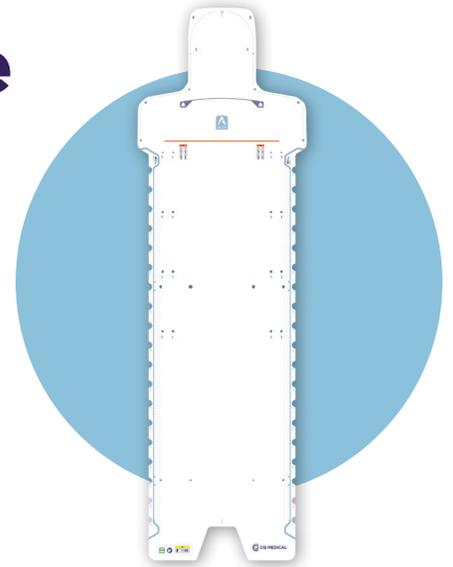

ALTA COMPATIBLE WITH HEXAPOD™ EVO RT SYSTEM

Enables full integration of CQ Medical positioning devices while maintaining clear access to the Elekta reference frame for attachment to the iBEAM® couchtop.



SEE ALTA
IN ACTION

Standard Package

INFINITYEDGE™ & INFINITYLOCK™

InfinityLock accessories tightly fasten to the device and can be infinitely indexed along the InfinityEdge of the device, which features Varian Style indexing.



ALTA ARMSHUTTLE™ ELITE

Quickly and securely attaches to provide comfortable, reproducible positioning of the arms out of the treatment field.



Alta Silverman Adapter Plate



Alta ArmShuttle Elite



InfinityLock Hand Grips




MR Safe Locating Bars



FlexLock Knee Support



FlexLock Foot Support

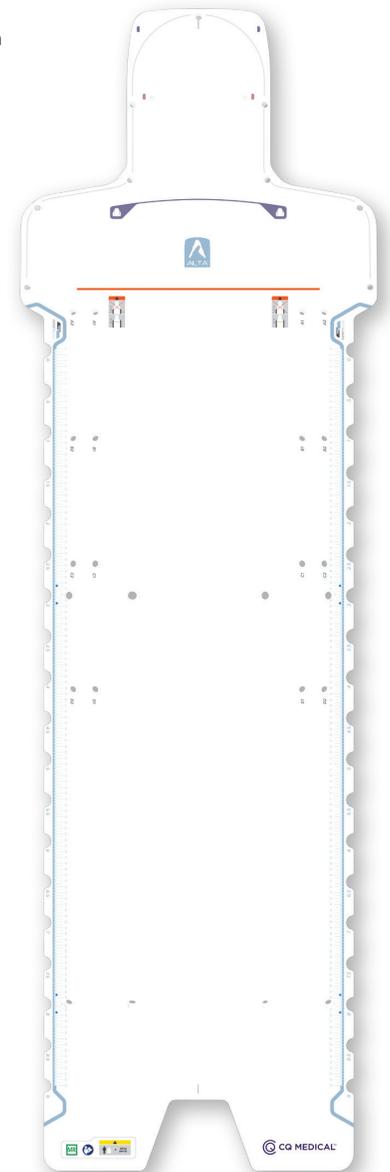

Wall Mount

Alta Standard Package

Alta Multipurpose Device	Qty 1
Alta ArmShuttle Elite	1
Alta Silverman Adapter Plate	1
FlexLock Knee Support	1
FlexLock Foot Support	1
InfinityLock Hand Grips	2
MR Safe Locating Bars	2
Wall Mount	1
	RT-4800-100

RT-4800-05 Alta Multipurpose Device for HexaPOD evo RT System

MTWBLO07 Alta Indexing Bar



Attach Two-Pin Wing Boards™ with the Alta Indexing Bar.